


## Soft Ride Measuring

- Properly measuring and fitting your horse for Soft Ride boots ensures they have maximum therapeutic effect and reduces the risk of the boots falling off or causing rubs.
- It is best to measure the hooves shortly after they have been trimmed, this will help ensure the boot is not too big. Regular hoof trimming every 5-6 weeks is ideal for horses with hoof issues, with a maximum of 8 weeks.
- Measure each foot individually, as some horses require two different sizes. Some horses may fall in between sizes.
- Please consult the attending veterinarian for sizing recommendations after measuring.
- After receiving the boots- hold the gel orthotic up to the bottom of the hoof. It should overlap the hoof by up to 1/8<sup>th</sup> on either side. If the foot overlaps the gel in length or width, then the boot is too small.

**Step 1: Pick out and clean the hooves intended for measuring**

**Step 2: Place the hoof squarely on a sturdy piece of paper or cardboard**

**Step 3: Using a pencil or pen, trace each hoof individually. Trace around the back of the hoof including the whole ground-bearing surface. Keep tracing close to the hoof wall or shoe.**

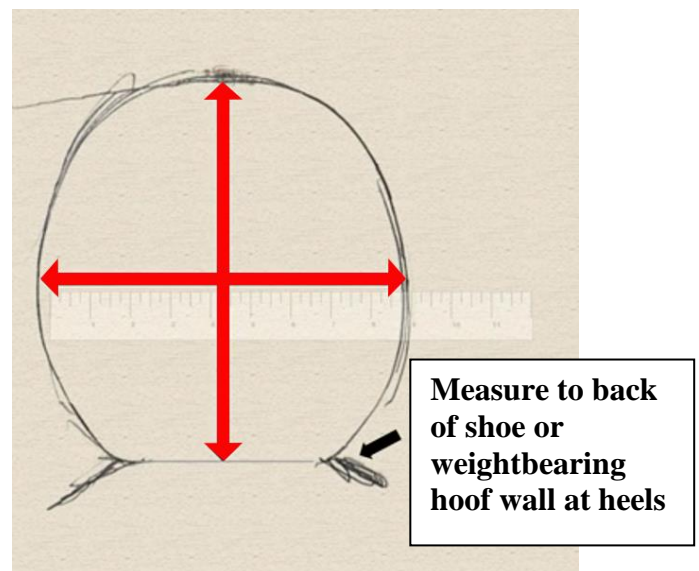

10995 St Rt 128 Harrison, OH  
 Phone # 513-715-2600  
 info@miamitownequine.com  
 miamitownequine.com



Anna Hood, DVM, CVMMP  
 Jamie Nickels, DVM

<b>SIZE</b>	<b>WIDTH Inches / mm</b>	<b>LENGTH mm</b>
# 2 Standard	6 5/8 in / 168mm	6 1/2 in / 164mm
# 2 Narrow	5 3/4 in / 146mm	6 3/8 in / 162mm
# 3 Long Narrow	5 3/4 in / 146mm	6 3/4 in / 171mm
# 3 Standard	6 3/8 in / 162mm	6 in / 152mm
# 3 Narrow	5 1/4 in / 133mm	6 in / 152mm
# 4 Standard	5 3/4 in / 146mm	5 in / 127mm
# 5 Long Standard	5 3/4 in / 146mm	5 1/2 in / 140mm
# 5 Long Narrow	4 3/4 in / 121mm	5 3/4 in / 146mm
# 5 Standard	5 3/4 in / 146mm	5 1/4 in / 133mm
# 6 Long Standard	5 3/8 in / 136mm	5 1/4 in / 133mm
# 6 Long Narrow	4 3/4 in / 121mm	5 1/4 in / 133mm
# 6 Standard	5 in / 127mm	4 1/2 in / 114mm
# 6 Narrow	4 1/4 in / 108mm	4 3/4 in / 121mm
# 7 Long Standard	4 7/8 in / 125mm	5 in / 127mm
# 7 Long Narrow	4 3/8 in / 110mm	5 in / 127mm
# 7 Standard	4 3/8 in / 110mm	4 in / 102mm
# 7 Narrow	4 in / 102mm	4 in / 102mm
# 8 Standard	3 3/4 in / 96mm	3 3/4 in / 96mm
# 9 Standard	3 1/4 in / 83mm	3 1/4 in / 83mm
# 10 Standard	2 3/4 in / 70mm	2 1/2 in / 64mm